

NEW


Advazorb

Areola



The protective and secure breast dressing developed with clinicians

- No more 'patchwork' dressings
- Low friction backing to prevent uncomfortable snagging on clothes
- Reduced dressing time and associated patient discomfort
- Central hole to prevent pressure on the nipple including incisions to accommodate larger nipples
- Designed for patient comfort
- Silfix soft silicone adhesive for atraumatic removal

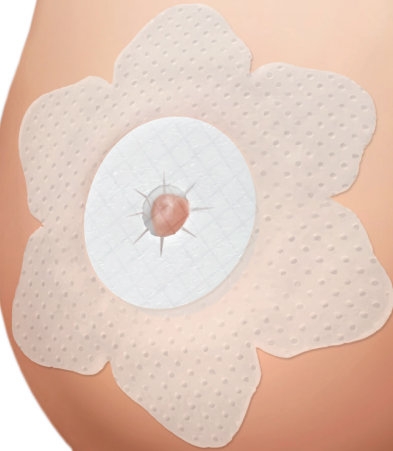


An Innovative Breast Dressing

Advazorb Areola is the result of a collaboration with Nottingham University Hospitals and Nottingham Breast Institute to fulfil a very specific need – a protective and secure breast dressing.

Advazorb Areola combines Silfix® soft silicone with conformable, absorbent foam to deliver excellent patient comfort, effective exudate management, secure fixation and atraumatic removal.

Three layers of performance

- 1 Low friction backing**
Increases 'wearability' and prevents uncomfortable snagging on clothes
- 2 Silfix silicone adhesive**
For easy placement and replacement
Atraumatic removal
- 3 Comfortable, absorbent foam**
With nipple hole to prevent nipple pressure.
Incisions ensure the hole suits larger nipple sizes.

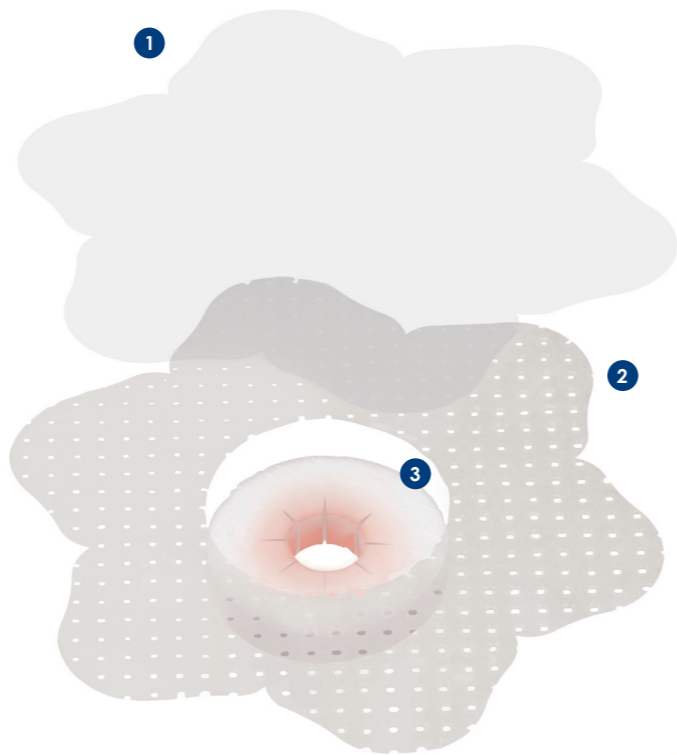


The challenge of dressing the breast

Clinicians are often faced with trying to create their own breast dressing from a patchwork of square dressings that do not conform to the shape of the breast.

Fit for purpose

Advazorb Areola is ergonomically designed to fit the contours of the breast making dressing application, easy, quick and comfortable for the patient.

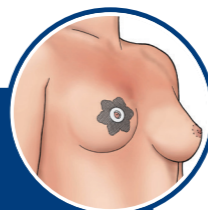

The ideal choice for:

- Mastopexy
- Breast reduction
- Nipple reconstruction
- Breast augment through areola incision
- Gynaecomastia
- Other Nipple / Areola surgery



Dimensions:

Total width **12.7cm**

Foam diameter **5.0cm**

Nipple hole **1.2cm plus incisions**

Product Code:

CR4490

For further information or to request a demonstration on the innovative Advazorb Areola, contact Advancis Medical on **+44 (0)1623 751 500** or email **Info@advancis.co.uk**



Find us online



www.advancis.co.uk



MAR605/V1